

www.saintyco.com

nspection Machine





CATALOGUE

Inspection Machine


The capsule flows from the hopper into a vibratory chute where a perforated screen eliminates coarse dust and chips, which are collected in a container. The capsules enter the inspection stage in a single layer. The capsules from the vibratory chute are placed on rotating rollers and are continuously rotate & transfer in front of the operator, which ensures a reliable all-sided visual inspection. An adjustable mirror enables the operator to inspect the opposite sides of the products. The rotation speed of the inspection rollers is adjustable individually, independent of the inspection speed.




I-500


Features:


- The rotation and inspection speed of capsule is infinitely adjustable. No hidden corners and blind areas in inspection stage.
- Inspecting rollers are made by FDA accepted transparent high density polymers. Lighting fixture installed underneath the roller conveyor belt, for illuminating the product from below, which is more suitable for inspection of transparent hard and softgel capsules.
- Defective capsules can be eliminated by the operator with a flexible handheld vacuum nozzle.
- Foot pedal switch to stop the conveyor belt at any time. While foot pedal is pressed, the products keep rotating in front of the operator at a separately adjustable speed.
- Lighting source from top and installed underneath the roller conveyor belt for illuminating the inspection stage in non-flashed fluorescent lamp.
 Fatigue of the vision can be reduced.
- For cleaning purposes the machine can be completely dismantled without any tools within a few minutes, easy to clean, no cross contamination risk.
- No noisy external vacuum cleaner necessary through the use of compact Venturi suction system Reduced fatigue of the operator through ergonomic design.

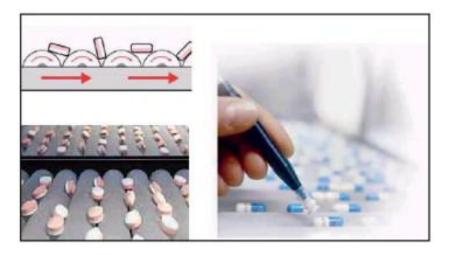

Dimensions:

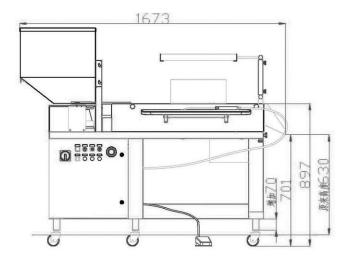

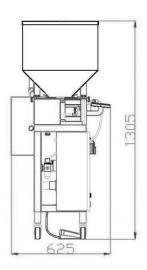

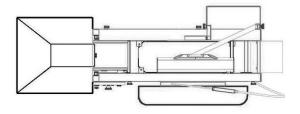

Technical Data	I-500
Application	Inspection of hard or soft capsules plain or coated tabletsand other
	round or cylindrical objects.
Min Tablet Thickness	3mm
Output	Up to 120, 000 capsules/h
Feed hopper	30 L
Electrical Connection	220V/50Hz, single phase
Power Input	0.6 kw
Compressed Air	0.6 MPa
Dimensions	1850 x 850 x 1315 mm
Weight	125 kg

SaintyCo also offer:

Encapsulation Solutions

Tablet Compression Solutions

Coating Systems

Cartoning Systems

Blister Packing Solutions

