



TRUSTAR PHARMA & PACKING EQUIPMENT CO., LTD.


TEL: 86-577-65905236 FAX: 86-577-65905150

Mobile: 86-13958877535

E-mail: tracy@trustar-pharma.com

[Http://www.trustar-pharma.com](http://www.trustar-pharma.com)

YH25 Series Powder Forming Machine

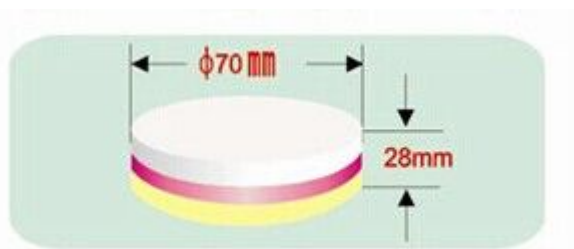

1. Machine Description

This machine is mainly used in pharmaceutical industry and also the chemical, food, electronic, chemical fertilizer and ceramic industries.

This machine can press the various kinds of granular materials or powder materials into tablets, such as round shape tablets, irregular shape tablets, double side engraved tablets and ring type tablets. This machine also can press single layer tablets, double layers tablets or triple layers tablets to meet the requirement of different customers.

2. Some technical features



TRUSTAR PHARMA & PACKING EQUIPMENT CO., LTD.

TEL: 86-577-65905236 FAX: 86-577-65905150

Mobile: 86-13958877535

E-mail: tracy@trustar-pharma.com

[Http://www.trustar-pharma.com](http://www.trustar-pharma.com)

- 1) Made of stainless steel, the housing is fully closed. The surface of the rotary turret is covered with a hardened layer so that the turret surface is wear resistant. The machine conforms to GMP requirements.
- 2) Adopt transparent windows, tableting state can be observed clearly. The windows can be opened, cleaning and maintenance is easy.
- 3) Machine design is reasonable. Operation, dismantling and maintenance is convenient.
- 4) PLC and touch screen control, all the running parameter can be set and shown.
- 5) The machine can stop automatic in case of over pressure.
- 6) all the drive device are located inside the machine so keep the machine clean.

3. Main technical parameters

Number of Station	25
Max. Main Pressure	250KN
Max. Tablet Diameter	70mm
Max. Filling Depth	60mm
Max. Tablet Thickness	25mm
Max. Turret Speed	16 (r/min)
Max. Production Capacity	24000(pcs/h)
Motor Power	11kw
Overall size	1520*1480*2600mm
Machine Weight	8500kg

4. Components List

No.	Name	Specification	Qty.	Brand
1	switchgroup	HZ12-40 40A	1PCS	DELIXI
2	breaker	DZ47-63 3P40A	1PCS	DELIXI



TRUSTAR PHARMA & PACKING EQUIPMENT CO., LTD.

TEL: 86-577-65905236 FAX: 86-577-65905150

Mobile: 86-13958877535

E-mail: tracy@trustar-pharma.com

[Http://www.trustar-pharma.com](http://www.trustar-pharma.com)

3	breaker	DZ47-63 3P10A	1PCS	DELIXI
4	breaker	DZ47-63 1P10A	1PCS	DELIXI
5	breaker	DZ47-63 1P4A	2PCS	DELIXI
6	breaker	DZ47-63 1P3A	2PCS	DELIXI
7	midget relay	HH54P/TP514X DC24V	1PCS	DELIXI
8	contactor	CDX1-10/OBB4 DC24V	1PCS	Shanghai machine tool equipment factory
9	contactor	CJX1-16 DC24V 16A	1PCS	DELIXI
10	power switch	NES-150-24	1PCS	Shangyuan
11	transducer	MD280NT15GB/18.5PB 15KW/18.5KW	1SET	Inovance
12	touch panel	1T5070 24VDC	1PCS	Inovance
13	communication cable	One meter	1PCS	Inovance
14	PLC	H1U-0806	1PCS	Inovance
15	scram button	LAY8-11ZS/R	1PCS	DELIXI
16	axial flow fan	125F2Y2 25W	2PCS	DELIXI
17	proximity switch	J3-D4C1 DC24Vwith 5-meter Cable plug-in line、10-36VDC、200MA	6PCS	DELIXI
18	plug and socket	013/313	1SET	SUNIX
19	socket	314	3PCS	SUNIX
20	plug and socket	63A-6H/380V-415	1SET	TY
21	plug and socket	WS55J40TP	2SET	TY
22	plug and socket	WP20J4TP WP20K4Z	1SET	TY
23	Terminal Block	made based on drawing	1PCS	Wanke, Tianjin



TRUSTAR PHARMA & PACKING EQUIPMENT CO., LTD.

TEL: 86-577-65905236 FAX: 86-577-65905150

Mobile: 86-13958877535

E-mail: tracy@trustar-pharma.com

[Http://www.trustar-pharma.com](http://www.trustar-pharma.com)

24	Install the soft copper wire	RV-32/0.2 1 square red 100 meters	meter	NH-YJV
25	Install the soft copper wire	RV-48/0.2 1.5 square black 100 meters	meter	NH-YJV
26	Install the soft copper wire	RV-48/0.2 1.5 square yellow-green 50 meters	meter	NH-YJV
27	Install the soft copper wire	RV-77/0.2 4 square black 50 meters	meter	NH-YJV