


RJWJ-200C Softgel Encapsulation Machine

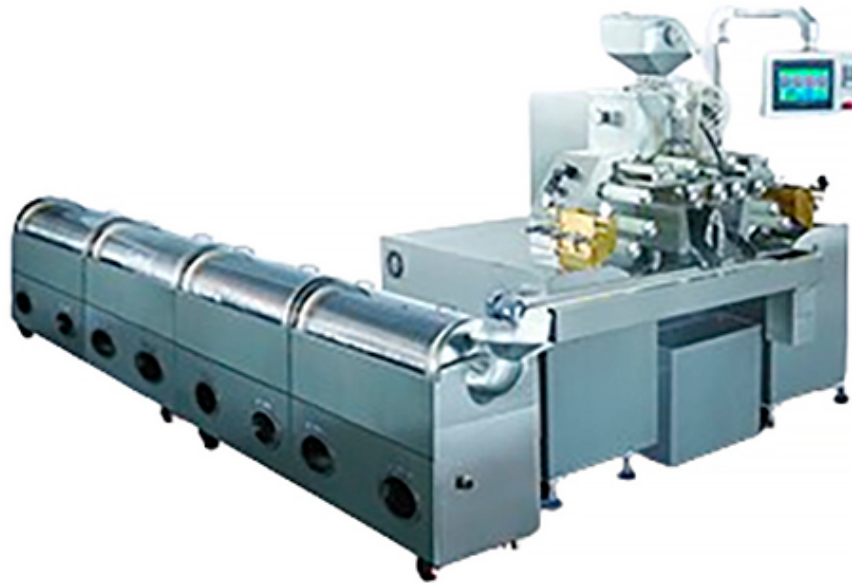

Softgel encapsulation machine can quantification ally inject the oil, mixing suspension and pasty mass into the capsule to form all kinds of softgel encapsulation with different size, shape and color. Because of the merits of high speed disintegration, easy to oral, hard to oxygenation, easy to store and take, the softgel encapsulation is welcomed by more and more customers. Owing to the high productivity, charging exactly, quality stabilization, high finished product rate, it is favored by more and more manufacturers and used in the fields of medicine, health products, cosmetic and paintball and so on.

This product has passed provincial-level appraise of new product and the certification of EU, adopted two items of patent technology in manicuring with exquisite workmanship.

Characteristic:

- ◆ Strong product capacity, it can produce 300 million granules (No.8 oval shape, 3 revolutions/minute, 24 hours / day, 300 days/year).
- ◆ Injecting medicine part adopts advanced linearity driving system, and it makes sure that the precision of charging and filling is less than 2%. Capsule package rate is up to 98%.
- ◆ The gelatin ribbon is lubricated by fine oil supply technology in order to lower the consumption of lubrication oil. The ribbon is lubricated by vegetable oil, so the capsules produced is free of wash.
- ◆ High quality capsule deferent equipment adopts food conveyor belt authorized by American FDA, without poison and adhibition and easy to clean.
- ◆ Optimized designed mould is made of high quality aviation alloy material by high precision numerical control process to insure the mould's service life.
- ◆ Advanced process and assembly technics make the complete appliance quality stabilization , high working rate, and good shape pill.
- ◆ Adopt automatic system project of supplying capsule and material to reduce workers' labor intensity and improve the working environment, which accords with the requirement of GMP.
- ◆ The liquid level in the gelatin box is fixed automatically in the constant position, and the centigrade mater display the thickness of rubber, which is easy and accurate in adjustment.

- ◆ The power of main motor is 1.5KW, which can avoid the phenomenon of insufficient driving power when the content is rather thick (take 185 model for example).
- ◆ Adopt the nonpoisonous materials to human body or now dropping matters into medicine liquid in this machine and medicine liquid, and the touching parts for gelatin.

Technical Parameters

Model	RJWJ-200C
Items	
Rotary speed of die roll	0.5rpm
Supplying volume for flow single-pillar	0-2ml
Content uniformity	±2%
Power for main machine	7Kw
Pillar quantity	12
Overall size of main machine	1500*950*2000mm
Weight of main machine	1000Kg